



NOVAmedtek

**NOVAMAG
NT-60**

**Magnetic Incontinence
Therapy Chair**


- The **NOVAmag NT-6o** is a luminous solution that is developed for treating urinary Incontinence by delivering pulsed magnetic waves to the weakened pelvic floor area. The pulsed magnetic waves depolarize the pelvic floor nerves and muscles, providing, in essence, pelvic muscle exercises without effort from the patient who is sitting on the treatment chair.
- The therapy carried out by **NT-6o** can act as a useful addition to conservative therapies as it is free of complications. The patient is able to receive the treatment fully clothed, since the she/he doesn't require prior preparation. The magnetic field passes through clothing and patient's tissue easily. The treatment is non-invasive, painless and with no side effects recorded against the magnetic therapy by **NT-6o**. Thus by feeling the magnetic waves, the patient can focus and feel the awareness of the muscle contractions to actively perform the exercises simple and in an effective way.
- A typical Incontinence treatment with **NT-6o** is 15 sessions, in which 3 sessions per week each for 20 minutes can be provided. After the first adaptations of the sessions, improvements will be observed, however, completion of the treatment course is must for advanced response under the consultant of the physician.

The NOVAmag NT-6o has also proven results in the following symptoms:

- Rehabilitation after pelvic floor surgery
- Stress, urge and mixed urinary incontinence
- Vaginal toning after delivery
- Dysmenorrhea (menstrual pain)
- Female orgasmic dysfunction
- Faecal incontinence
- Postoperative urinary incontinence
- Chronic prostatitis inflammation

SYSTEM FEATURES:

Therapy Frequency : 10 - 50Hz
 Operating Power : 30-100 %
 Operating Mode : High – Low Frequency Therapies
 Therapy Time : Adjustable 5-40 Minutes
 Magnetic Wave : Up to 3 Tesla

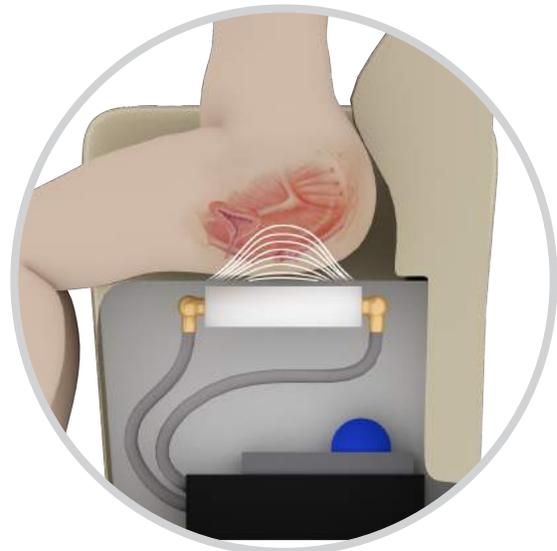

LEATHER COLOURS:



**Leather Colours can be chosen according to customer's request.*